

User Manual

Luna® II CPAP / Luna® II Auto CPAP

LG2000 / LG2A00



Table of Contents

1. Introduction	1
2. Symbols	1
2.1 Control Buttons	1
2.2 Device Symbols	1
3. Warning, Caution and Important Tip	4
4. Intended Use	4
5. Contraindications	5
6. Specifications	6
7. Available Therapies	8
8. Glossary	8
9. Model	9
10. Package Contents	10
11. System Features	11
12. First Time Setup	12
12.1 Placing the Device	12
12.2 Installing the Air Filter and Filter Cap	13
12.3 Connecting to Power	13
12.4 Assembling the Tube and Mask	14
12.5 Using Oxygen with the Device	15
12.6 Inserting the SD Card (Only for the device that equipped with SD card)	17
12.7 Using the Cellular Module	18
12.8 Starting Treatment	18
13. Routine Use	18
13.1 Connecting the Tube	18
13.2 Adjusting the Tube	18
13.3 Turning on the Airflow	18
13.4 Heating the Water	18
13.5 Using the Ramp Button	19
13.6 Turning the Device Off	19
14. Heated Humidifier	19
14.1 Filling the Water Chamber	19
14.1.1 Removing the Water Chamber	19
14.1.2 Filling Water	20
14.1.3 Inserting the Water Chamber	21
14.2 Emptying the Water Chamber	21
14.3 Setting the Humidity Level	22
15. Navigating the Patient Menu	23
15.1 Steps to Navigating the Patient Menu	23
15.1.1 Accessing the Main Interface	23
15.1.2 Bringing up the Initial Setup Interface	23
15.1.3 Accessing the Setup Interface	24
15.1.4 Selecting Options	24
15.1.5 Adjusting Options	24
15.1.6 Confirming Adjustments	25
15.1.7 Turning Pages	25
15.1.8 Exiting the Patient Menu	25

15.2 Options of the Patient Menu and Corresponding Descriptions	27
16. Alert	28
17. Cleaning	29
17.1 Cleaning the Mask and Headgear	30
17.2 Cleaning the Cellular Module	30
17.3 Cleaning the Water Chamber	30
17.4 Cleaning the Enclosure	31
17.5 Cleaning the Tube	31
17.6 Cleaning and Replacing the Air Filter	31
18. Traveling	32
18.1 Traveling with the Device	32
18.2 Traveling by airplane	32
19. Reordering	32
20. Technical Support	33
21. Disposal	33
22. Troubleshooting	33
22.1 Common Problems in Patients and Corresponding Solutions	34
22.2 Common Problems in the Device and Corresponding Solutions	36
23. EMC Requirements	37
24. Limited Warranty	43

1. Introduction










Thank you for your purchase of the Luna® II CPAP / Luna® II Auto CPAP. This User Manual will provide you with information about your device. Please read it carefully. If you experience any difficulties or problems during use, please contact your homecare provider or physician.

2. Symbols

2.1 Control Buttons

	Ramp Button
	Mute Button
	Knob

2.2 Device Symbols

	Follow Instructions for Use
	Consult Instructions for Use
	Type BF Applied Part (mask)
	Class II (Double Insulated)
	AC Power
	DC Power
IP22	≥ 12.5 mm Diameter, Dripping (15° tilted)
	Hot Surface
	Logo of BMC Medical Co., Ltd.
	Serial Number



Manufacturer



Disassembly is prohibited



Batch Code



Non-ionizing Radiation



SD Card



Marking of Electrical and Electronic Equipment



Air Inlet



Air Outlet



Made in China



MR Unsafe



Complies with RTCA DO-160 section 21, category M.



Caution



Medical Device



Prescription only



Catalogue number



Unique device identifier



Model Number



Temperature Limit



Humidity Limit



Atmospheric Pressure Limit

3. Warning, Caution and Important Tip

WARNING!

Indicate the possibility of injury to the user or operator.

CAUTION!

Indicate the possibility of damage to the device.

IMPORTANT TIP!

Indicate the possibility that such operation may affect the effectiveness or ease of use of the device.

Warnings, Cautions, and Important Tips appear throughout this manual as they apply.

4. Intended Use

The Luna® II CPAP and Luna® II Auto CPAP are intended to deliver positive pressure for the treatment of obstructive sleep apnea. The optional integrated heated humidifier is indicated for the humidification and warming of air from the flow generator. These devices are intended for single patient use by prescription in the home or hospital/institutional environment on adult patients.

WARNINGS!

- This device is intended for adult use only.
- This device is not intended for life support.
- The instructions in this manual are not intended to supersede established medical protocols.
- Do not bring the device or accessories into a Magnetic Resonance (MR) environment as it may cause unacceptable risk to the patient or damage to the device or MR medical devices. The device and accessories have not been evaluated for safety in an MR environment.
- Do not use the device or accessories in an environment with electromagnetic equipment such as CT scanners, Diathermy, RFID and electromagnetic security systems (metal detectors) as it may cause unacceptable risk to the patient or damage to the device. Some electromagnetic sources may not be apparent, if you notice any unexplained changes in the performance of this device, if it is making unusual or harsh sounds, disconnect the power cord and discontinue use. Contact your home care provider.

CAUTIONS!

- Federal law restricts this device to sale by or on the order of a physician.
- The device is intended for use by operators trained or experienced in similar equipment.
- The patient is an intended operator.

- Cleaning can be performed by the patient.

IMPORTANT TIP!

- Read and understand the entire user manual before operating this system. If you have any questions concerning the use of this system, contact your home care provider or health care professional.

5. Contraindications

Studies have shown that the following pre-existing conditions may contraindicate the use of positive airway pressure therapy for some patients:

Absolute Contraindications: Pneumothorax, mediastinal emphysema; cerebrospinal fluid leak, traumatic brain injury, or pneumocephalus; shock caused by a variety of conditions before treatment; active epistaxis; upper gastrointestinal bleeding before treatment; coma or impaired consciousness making the use of mask during therapy impossible; giant vocal fold polyp, etc.

Relative Contraindications: Severe coronary heart disease complicated with left ventricular failure, acute otitis media, excessive respiratory secretions and weak cough, weak spontaneous breathing, nasal or oral tracheal intubation and tracheotomy, severe nasal congestion caused by a variety of conditions, lung bullae, and allergies to breathing masks, etc.

The following side effects may occur during treatment:

- Dryness of the mouth, nose and throat
- Abdominal bloating
- Ear or sinus discomfort
- Eye irritation
- Skin irritation due to the use of a mask
- Chest discomfort

IMPORTANT TIPS!

- An irregular sleep schedule, alcohol consumption, obesity, sleeping pills, or sedatives may aggravate your symptoms.
- Please use a mask which meets ISO 17510:2015.

CAUTION!

- Contact your health care professional if symptoms of sleep apnea recur. Contact your health care professional if you have any questions concerning your therapy.

6. Specifications

Device Size

Dimensions: 274 mm × 184 mm × 115 mm

Weight: 1.9 kg

Water capacity: To maximum fill line 360 mL

Product Use, Transport and Storage

Operation	Transport and Storage
Temperature: 5°C to 35°C (41°F to 95°F)	-25°C to 70°C (-13°F to 158°F)
Humidity: 15% to 93% Non-condensing	15% to 93% Non-condensing
Atmospheric Pressure: 760 to 1060 cmH ₂ O	760 to 1060 cmH ₂ O
Altitude: Sea level to 2300 m	Sea level to 2300 m

Heated Humidifier

Humidifier Settings: off, 1 to 5 (95°F to 154.4°F / 35°C to 68°C)

Humidifier Output: No less than 10 mg H₂O/L

Environmental Conditions: Maximum airflow, 35°C, 15% relative humidity

Mode of Operation

Continuous

Work Mode

For LG2000: CPAP

For LG2A00: CPAP, Auto

SD Card

The SD card can record patient data and fault information

AC Power Consumption

100 V - 240 V \sim , 50 Hz / 60 Hz, Max 2 A

Device offer to Communications Port

5 V $\overline{\text{---}}$ 2.0 A

Type of Protection Against Electric Shock

Class II Equipment

Degree of Protection Against Electric Shock

Type BF Applied Part

Degree of Protection Against Ingress of Water

IP22

Pressure Range

4 to 20 cmH₂O (in 0.5 cmH₂O increments), ≤ 30 cmH₂O under single fault conditions.

Static Pressure Stability

±0.5 cmH₂O

Ramp

The ramp time ranges from 0 to 60 minutes.

The A-weighted sound pressure level and sound power level

When the device is working at the pressure of 10 hPa, its sound pressure level and sound power level shall not be greater than the values in the following table:

Sound Pressure Level	Uncertainty	Sound Power Level	Uncertainty
28 dB(A)	2 dB(A)	36 dB(A)	2 dB(A)

Note: Declared dual-number noise emission values in accordance with ISO 4871:1996.

Maximum Flow

Test Pressures (cmH ₂ O)	4	8	12	16	20
Measured Pressure at the Patient Connection Port (cmH ₂ O)	3	7	11	15	19
Average Flow at the Patient Connection Port (L/min)	90	150	150	150	145

When the working pressure is set to the values listed in the table, the average flow rate at the patient end should be greater than 80% of the corresponding flow value in the table.

Automatic leak compensation

The pressure error $\leq \pm 1$ cmH₂O when the air leakage is no more than 70 L / min.

Pressure

Range: 0 ~ 20 cmH₂O

Margin of Error: $\pm (0.4 \text{ cmH}_2\text{O} + 4\%)$

Tube

Length: 6 ft. (1.8 m)

Maximum Delivered Gas Temperature


$\leq 43^\circ\text{C}$

The Form and the Dimensions of the Patient Connection Port

The 22 mm conical air outlet complies with ISO 5356-1.

7. Available Therapies

The device delivers the following therapies:

CPAP – Delivers Continuous Positive Airway Pressure; CPAP maintains a constant level of pressure throughout the breathing cycle. If your health care professional has prescribed ramp for you, you can press **the Ramp Button**  to reduce the pressure and then gradually increase the pressure to the therapeutic pressure setting so that you can fall asleep more comfortably.

Auto – Delivers CPAP therapy and provides an air pressure no less than the prescribed one based on the patient's needs.

8. Glossary

Apnea

A condition marked by the cessation of spontaneous breathing.

Auto

Adjust CPAP pressure automatically to improve patient comfort based on monitoring of apnea and snoring events.

Auto Off

When this feature is enabled, the device automatically discontinues therapy whenever the mask is removed.

Auto On

With this feature, the device automatically initiates therapy when you breathe into the mask.

CPAP

Continuous Positive Airway Pressure.

iCode

A feature that is intended to give access to compliance and therapy management information. The "iCode" consists of six separate codes displayed in the Patient Menu, each code is a sequence of numbers. The "iCode QR" and "iCode QR+" display two-dimensional codes.

LPM

Liters Per Minute.

OSA

Obstructive Sleep Apnea.

Patient Menu

The display mode in which you can change patient-adjustable device settings, such as the starting pressure for the Ramp feature.

Ramp

A feature that may increase patient comfort when therapy is started. It begins from a low pressure and then gradually increases to the prescribed setting pressure so that the patient can fall asleep more comfortably.

Reslex

A therapy feature that is enabled by your home care provider to provide pressure relief during exhalation.

Standby State

The state of the device when power is applied but the airflow is turned off.

min

Means the time unit "minute".

h

Means the time unit "hour".

yy mm dd / mm dd yy / dd mm yy

Denotes date.

9. Model

Item No.	Product Description			
	Product Contents	Optional Accessory	Work Mode	Maximum Work Pressure (cmH ₂ O)
LG2000	Device (2.4-inch TFT)	Tube, Mask, Cellular Module,	CPAP	20
LG2A00	Device (3.5-inch TFT)		CPAP, Auto	

10. Package Contents

After unpacking the system, make sure you have everything shown here (Different models of the product may contain different components):

No.	Articles	Qty.	Notes
1	Device	1	
2	Air Filter	2	
3	Power Adapter	1	
4	Power Cord	1	
5	Tube	1	
6	SD Card	1	Optional
7	Carrying Case	1	
8	Enclosed Documents	1	

All parts and accessories are not made with natural rubber latex.

If the use, maintenance, cleaning and disinfection are in strict accordance with the User Manual:

The expected service life of the device is 5 years from first date of use.

The expected service life of the Air Filter is 6 months from first date of use.

The expected service life of the Tube is 6 months from first date of use.

The expected service life of the Water Chamber is 6 months from first date of use.

IMPORTANT TIPS!

- If any of the above parts are missing, contact your home care provider.
- Contact your home care provider for additional information on the available accessories of this device. When using optional accessories, always follow the instructions enclosed with the accessories.

WARNINGS!

- This device should only be used with the mask and accessories manufactured or recommended by REACT HEALTH or with those recommended by your prescribing physician. The use of inappropriate masks and accessories may affect the performance of the device and impair the effectiveness of therapy.
- Do not pile up the long tubing at the head of the bed, as it may wrap around the head or neck of the patient during sleep.
- Do not connect any equipment to the device unless recommended by REACT HEALTH or your health care provider.
- Please contact REACT HEALTH to obtain an SD card if needed.
- Exceeding the Expected Service life, our company cannot guarantee the normal function of the device, nor the safety and effectiveness of the device.

11. System Features

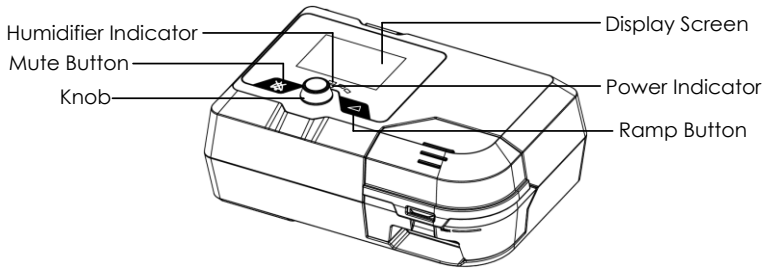


Fig. 11-1

Name	Function
Humidifier Indicator	Indicate the humidity level. There are five levels in total. The number of indicator lights that light up is directly proportional to the humidity level. If the indicator lights are off, it means the humidifier is turned off
Mute Button	Press this button to mute the alert. However, if the problem causing the alert is not solved, the alert will sound again two minutes later
Knob	Start treatment and adjust device settings
Ramp Button	Enable the Ramp feature
Display Screen	Display menus for operation, messages, monitoring data, etc.
Power Indicator	Indicate the power supply status

CAUTION!

- The pictures in this manual are only for reference, if they are different from the material object, the latter shall prevail.

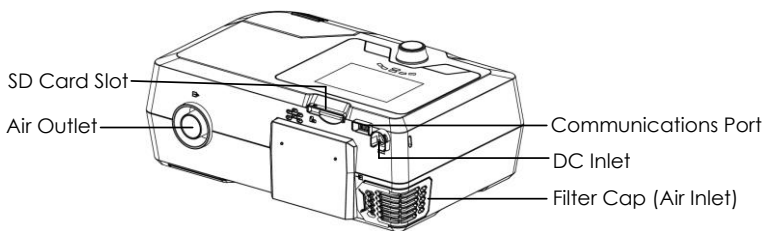


Fig. 11-2

Name	Function
SD Card Slot	Insert the SD card into this slot
Air Outlet	Deliver pressurized air; connects to the tube
Communications Port	Connected to external equipment (Not for connection to un-recommended devices)
DC Inlet	An inlet for the DC power supply
Filter Cap (Air Inlet)	Place the cap on the air filter, which is used to filter dust and pollen in the air entering the device

12. First Time Setup

12.1 Placing the Device

Place the device on a firm, flat surface.

WARNINGS!

- If the device has been dropped or mishandled, if the enclosure is broken, or if water has entered the enclosure, disconnect the power cord and discontinue use. Contact your home care provider immediately.
- If the room temperature is warmer than 95°F (35°C), the airflow produced by the device may exceed 109.4°F (43°C). The room temperature must be kept below 95°F (35°C) while the patient uses the device.

CAUTIONS!

- If the device has been exposed to either very hot or very cold temperatures, allow it to adjust to room temperature (approximately 2 hours) before beginning setup.
- Make sure the device is away from any heating or cooling equipment (e.g., forced air vents, radiators, air conditioners).
- The device is not suitable for use in high humidity environments. Make sure that no water enters the device.
- Make sure that bedding, curtains, or other items are not blocking the filter or vents of the device.
- Keep pets or children away from the device and avoid small objects being

inhaled or swallowed.

- To avoid explosion, this device must not be used in the presence of flammable gases (e.g. anesthetics).
- Tobacco smoke may cause tar build-up within the device, leading to the malfunctioning of the device.
- Air must flow freely around the device for it to work properly.
- Empty the water chamber completely before moving the device.

12.2 Installing the Air Filter and Filter Cap

(1) Attach the air filter to the filter cap, as shown in Fig. 12-1.

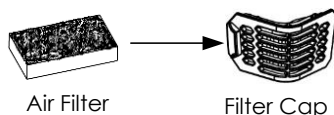


Fig. 12-1

(2) Install the filter cap containing the air filter to the device, as shown in Fig. 12-2.

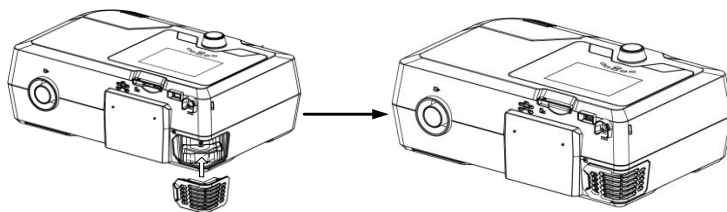


Fig. 12-2

CAUTIONS!

- The air filter must be in place when the device is operating.
- Installing the air filter and filter cap, device must be unplugged.
- Please change the Air Filter regularly and don't block it; Fire, open ignition source and smoking prohibited. (Refer to 17.6 Replacing the Air Filter)

12.3 Connecting to Power

- (1) Insert the plug of the power adapter into the DC Inlet on the back of the device;
- (2) Connect the power cord to the power adapter;
- (3) Plug the other end of the power cord into the power outlet.

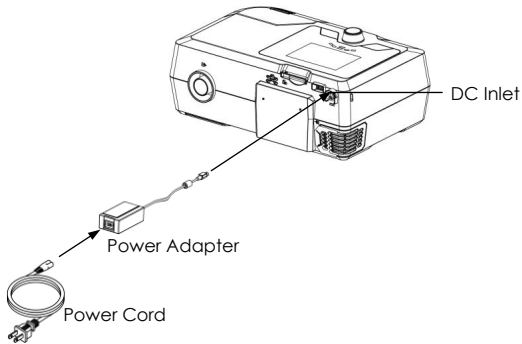



Fig. 12-3

WARNINGS!

- The device is powered on for use when the power cord and power adapter is connected. The **Knob**  turns the blower On / Off.
- Use of the device at an AC voltage beyond the stated range (see Section 5 "AC Power Consumption") may damage the device or cause device failure.
- Connect to appropriate power for proper operation of the device.

CAUTION!

- Inspect the power cord often for any signs of damage. Replace a damaged cord immediately.

IMPORTANT TIPS!

- After interruption and restoration of the power supply, the device will restore its pre-interruption working status automatically.
- To remove AC power, disconnect the power cord from the power outlet.

12.4 Assembling the Tube and Mask

(1) Connect one end of the tube to the air outlet of the device, as shown in Fig. 12-4.

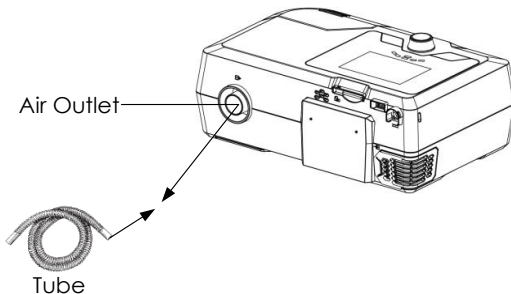


Fig. 12-4

(2) Connect the other end of the tube to the mask according to the user manual for the mask.

WARNINGS!

- If you are using a mask with a built-in exhalation port, connect the mask's connector to the tube.
- If you are using a mask with a separate exhalation port, connect the tube to the exhalation port. Position the exhalation port so that the vented air is blowing away from your face. Connect the mask's connector to the exhalation port.
- If you are using a full-face mask (a mask covering both your mouth and nose), the mask must be equipped with a safety (entrainment) valve.
- In order to minimize the risk of CO₂ rebreathing, the patient should observe the following instructions:
 - Use the accompanying tube and mask provided by REACT HEALTH.
 - Do not wear the mask for more than a few minutes while the device is not operating.
 - Use only masks with vent holes. Do not block or try to seal the vent holes in the exhalation port.
- If condensation appears in the tube, remove then drain the tube. Turn down the humidification.

12.5 Using Oxygen with the Device

Oxygen may be added at the mask connection. Please observe the instructions listed below when using oxygen with the device.

WARNINGS!

- Connect the oxygen tube to the oxygen inlet of the mask.
- The oxygen supply must comply with the local regulations for medical oxygen.
- Turn on the device before turning on the oxygen. Turn off the oxygen before turning off the device. Explanation of Warning: When the device is turned off, but the oxygen flow still exists, oxygen may accumulate within the device's enclosure and pose a fire hazard. Turning off the oxygen before turning off the device will prevent oxygen accumulation in the device and reduce the risk of fire. This warning applies to CPAP devices.
- Oxygen supports combustion. Keep the device and the oxygen container away from heat, open flames, any oily substances, or other sources of ignition. DO NOT smoke in the area near your Luna® II or the oxygen container.
- Sources of oxygen should be located more than 1 m from the device.
- When using oxygen with this system, a Pressure Valve must be placed in-line with the patient circuit between the device and the oxygen source. The pressure valve helps prevent the backflow of oxygen from the patient circuit into the device when the unit is off. Failure to use the pressure valve could result in a fire hazard.

- When the pressure valve is installed, the Auto On function on the device will not work. To turn on the device either breath into your mask or select the On/Off button.
- Supplemental oxygen must not be used while smoking or in the presence of an open flame.
- When using the device with an oxygen supply, check the following:
Starting therapy - ensure the device is on and blowing air before the oxygen supply is turned on.
Stopping therapy - ensure the oxygen supply is turned off first, then the device.
This will ensure oxygen does not accumulate within the device and create a risk of fire.
- Do not connect the device to an unregulated or high pressure oxygen source. The pressure of oxygen source does not exceed the work pressure of the device.

12.6 Inserting the SD Card (Only for the device that equipped with SD card)

Insert the SD card into the SD Card Slot, as shown in Fig. 12-5.

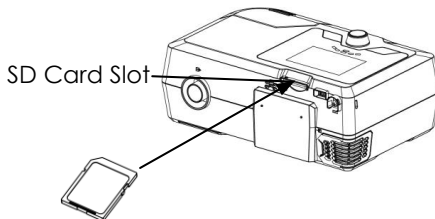


Fig. 12-5

If the SD card is inserted correctly, a symbol indicating correct insertion will appear in the Main Interface on the screen of the device, as shown in Fig. 12-6.

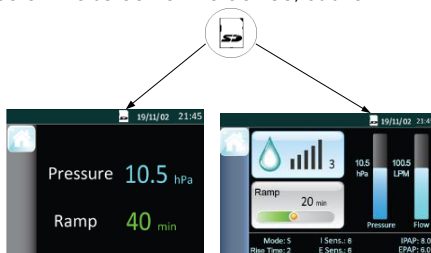


Fig. 12-6

If the SD card is inserted incorrectly, a symbol indicating incorrect insertion will appear in the Main Interface on the screen of the device, as shown in Fig. 12-7.

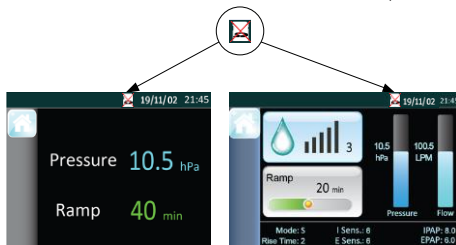


Fig. 12-7


CAUTION!

- If the SD card is not inserted, there will not be a symbol in the Main Interface on the screen of the device.
- To avoid data loss or any damage to the SD card, the SD card can only be removed after the device stops delivering air.

12.7 Using the Cellular Module

For more details, please refer to the corresponding user manual.

12.8 Starting Treatment

Connect the device to a power outlet, press **the Knob** , and the device will start delivering air.

WARNINGS!

- Be sure to follow your physician's instructions on adjusting the settings! To order any accessories not included with this device, contact your home care provider.
- DO NOT connect any ancillary equipment to this device unless recommended by REACT HEALTH or your physician. If you suffer from chest discomfort, shortness of breath, stomach bloating, or severe headache when using the device, contract your physician or qualified medical personnel immediately.

13. Routine Use

13.1 Connecting the Tube

Connect the power cord, power adapter, and tube properly according to the instructions in the First Time Setup (Chapter 12). Connect the mask and headgear according to the user manual for the mask.


CAUTION!

- Before each use, examine the tube for any damage or debris. If necessary, clean the tube to remove the debris. Replace any damaged tube. Make sure that the mask does not leak.

13.2 Adjusting the Tube

Lie down on your bed, and adjust the tube so it is free to move if you turn during sleep. Adjust the mask and headgear until you have a comfortable fit and until there are no airflow leaks around the mask.

13.3 Turning on the Airflow

Press **the Knob**  to turn on the airflow. The screen will display treatment pressure and other information.


13.4 Heating the Water

Pay attention to the humidifier indicator lights when using the humidifier. The indicator lights indicate the On / Off state of the humidifier. It is off when all indicator lights go out.

CAUTION!

- Observe the water level of the water chamber before using the humidifier. Make sure there is sufficient water in the water chamber, and avoid heating the device with an empty water chamber.


13.5 Using the Ramp Button

Every time **the Ramp Button**  is pressed, the pressure will drop to the initial pressure, and then gradually rise to the prescribed treatment pressure according to the preset ramp time, so as to make the patient fall asleep easily. The screen displays a real-time countdown of the remaining ramp time in minutes.

CAUTIONS!

- You can press **the Ramp Button**  as often as you wish during sleep.
- The ramp feature is not prescribed for all users.

13.6 Turning the Device Off

Take off the mask and headgear, press and hold **the Knob**  for two seconds, and the device will stop delivering air. Disconnect the power cord from the power outlet to power off the device.

CAUTIONS!

- Do not position the device where it is difficult to disconnect the power cord from the power outlet to power off the device.
- To isolate the device from the supply mains, disconnect the plug.

14. Heated Humidifier

The humidifier is available from your home care provider. The humidifier may reduce nasal dryness and irritation by adding moisture (and heat if applicable) to the airflow.

14.1 Filling the Water Chamber

14.1.1 Removing the Water Chamber

Grab the water chamber, and pull it out of the device, as shown in Fig. 14-1.

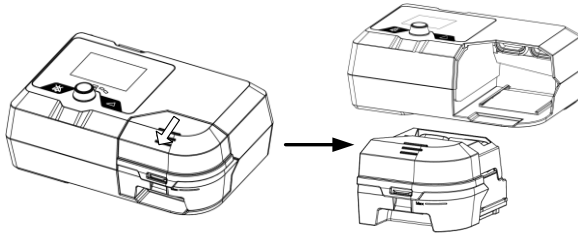


Fig. 14-1

WARNING!

- Turn the device off and allow approximately 15 minutes for the heater plate and water to cool.

14.1.2 Filling Water

Open the cap, as shown in Fig. 14-2, and fill the water chamber with approximately 360 mL of water, as shown in Fig. 14-3. Make sure that the water does not exceed the maximum water level line.

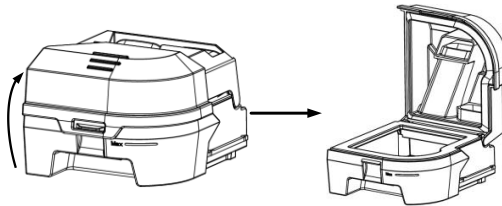


Fig. 14-2

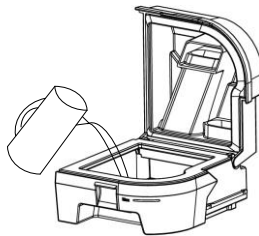


Fig. 14-3

WARNING!

- Change water before every use and do not surpass the MAX fill line.

CAUTIONS!

- Empty the water chamber when the heated humidifier is not in use.
- Distilled water is recommended.

14.1.3 Inserting the Water Chamber

Close the cap after it is filled with water, as shown in Fig. 14-4, and return it to the device, as shown in Fig. 14-5.

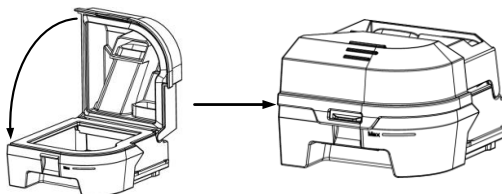


Fig. 14-4

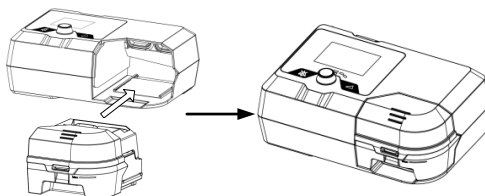


Fig. 14-5

WARNING!

- For safety purposes, the device must be placed on a flat surface at a level lower than the patient's head on a bed, so that the condensation flows back to the water chamber rather than remain in the tubing causing rainout.

CAUTIONS!

- Avoid moving or tilting the device when the water chamber has water in it.
- Take precautions to protect furniture from water damage.

14.2 Emptying the Water Chamber

(1) **Removing the water chamber** according to instructions in 14.1.1.

(2) **Emptying the water chamber:** Open the cap, as shown below, and pour any remaining water out of the water chamber.

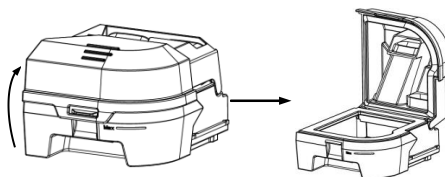



Fig. 14-6

CAUTION!

- Empty and air-dry the water chamber when the device is not in use.

(3) **Inserting the Water Chamber** according to instructions in 14.1.3.

14.3 Setting the Humidity Level

After the device is powered on, turn **the Knob**  to turn on or turn off the heated humidifier and to adjust the humidity level according to instructions on the screen of the device.

There are five humidity levels available, and the number of blue indicator lights that light up is directly proportional to the humidity level. If none of the indicator lights light up, it means that the heated humidifier is turned off.

The temperature of the water in the water chamber maintains a constant set level. Three indicator lights light up when the humidity is adjusted to Level 3, as shown in Fig. 14-7.

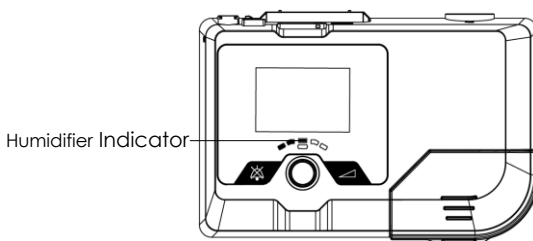


Fig. 14-7

CAUTIONS!

- Generally speaking, the humidity inside the mask is low when the water temperature is low.
- The greater the difference between the temperature inside the air tubing and room temperature is, the more easily condensation occurs inside the tubing.
- If there are only a few condensed water droplets inside the tubing in the morning after therapy, it means that the humidity level is appropriate; if there is lots of condensed water droplets inside the tubing and / or mask, it means that the humidity level is too high and should be set lower; Nasal dryness means that the humidity level is too low and should be set higher.

WARNING!

- Do not touch the heater plate of the device when it is working, otherwise you may get burned. Turn off the heat when the heated humidifier is not in use.

15. Navigating the Patient Menu

15.1 Steps to Navigating the Patient Menu

15.1.1 Accessing the Main Interface

Connect the power cord and power adapter properly. The screen displays the Main Interface shown in Fig. 15-1 (applies to LG2000), or the Main Interface shown in Fig. 15-2 (applies to LG2A00).



Fig. 15-1

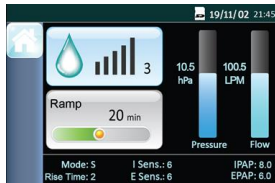



Fig. 15-2

15.1.2 Bringing up the Initial Setup Interface

From the Main Interface shown in Fig. 15-1 or Fig. 15-2, press and hold the **Ramp Button**  for three seconds. The screen displays the Initial Setup Interface of the Patient Menu, as shown in Fig. 15-3.

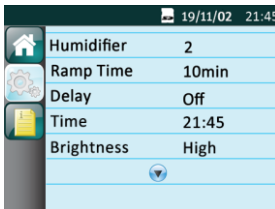






Fig. 15-3

The first icon  on the left side of the screen indicates the Main Interface, the second icon  indicates the Initial Setup Interface, and the third icon  indicates the iCode Interface. As you turn the **Knob** , the cursor switches among the three icons, and the interface displayed on the screen changes accordingly.

15.1.3 Accessing the Setup Interface





When the cursor is on the icon , the screen displays the Setup Interface. Access the Setup Interface by pressing **the Knob** . The first option on the Setup Interface is then displayed in blue, as shown in Fig. 15-4.



Fig. 15-4

15.1.4 Selecting Options

As you turn **the Knob**  clockwise, the cursor moves downwards from one option to another. As you turn it counterclockwise, the cursor moves upwards. When the cursor is on a certain option, press **the Knob** , and the option is then displayed in yellow, meaning that the option can now be adjusted, as shown by the **Humidifier** option in Fig. 15-5.

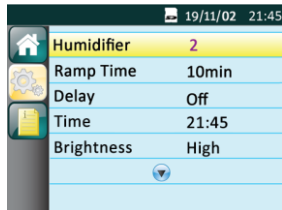





Fig. 15-5

15.1.5 Adjusting Options

Adjust the option by turning **the Knob** . As shown in Fig. 15-5, the **Humidifier** option is selected. As you turn **the Knob**  clockwise, the numbering increases, indicating a higher humidity level. As you turn **the Knob**  counterclockwise, the numbering decreases, indicating a lower humidity level. The **Humidifier** option is still displayed in yellow, as shown in Fig. 15-6.

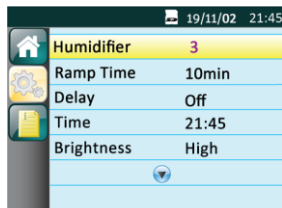



Fig. 15-6

15.1.6 Confirming Adjustments

Confirm your adjustment to an option by pressing **the Knob** . The option is then displayed in blue, as shown in Fig. 15-7.

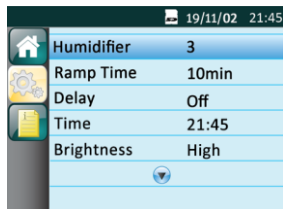



Fig. 15-7

15.1.7 Turning Pages

When the cursor is on **Brightness**, the last option shown in Fig. 15-7, the remaining options will appear on a new page if you continue to turn **the Knob**  clockwise, as shown in Fig. 15-8.

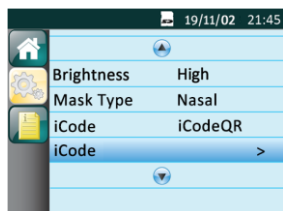


Fig. 15-8

Note:   are page turning symbols.

15.1.8 Exiting the Patient Menu

(1) Returning to the Initial Setup Interface

Move the cursor to the **Back** option by turning **the Knob** , as shown in Fig. 15-9.

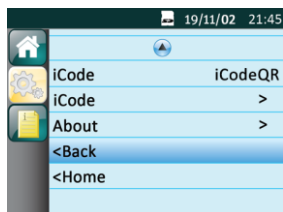




Fig. 15-9

Press **the Knob** , the cursor jumps to the second icon  on the left side of the screen. The screen displays the Initial Setup Interface, as shown in Fig. 15-10.

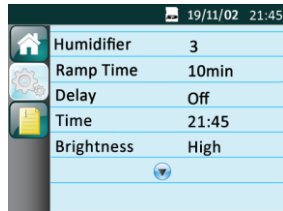



Fig. 15-10

(2) Returning to the Main Interface

Move the cursor to the **Home** option by turning **the Knob** , as shown in Fig. 15-11.

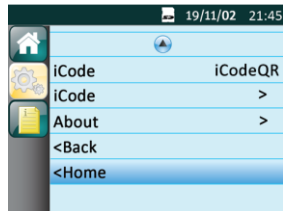


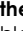



Fig. 15-11

Press **the Knob**  to exit the Patient Menu. The screen will display the Main Interface shown in Fig. 15-1 or Fig. 15-2.

15.2 Options of the Patient Menu and Corresponding Descriptions

Option	Range	Description
Humidifier	Off, 1 ~ 5	There are five humidity levels available. As the numbering increases, the humidity rises accordingly. "Off" means the humidifier is turned off.
Reslex	Off, 1 ~ 3	This feature enables the device to automatically reduce the treatment pressure when the patient exhales, so as to make the user more comfortable. The higher the numbering is, the more pressure the device reduces. "Off" means this feature is disabled.
Ramp Time	0 - Max Ramp	In order to increase comfort and help the patient fall asleep easily, the pressure can increase gradually, when the Ramp feature is enabled. The ramp time during which the initial pressure rises to the prescribed treatment pressure can be adjusted. As you turn the Knob  to the nearest point, the numbering increases or decreases by five minutes. The screen displays a real-time countdown of the remaining ramp time in minutes
Delay	On / Off	When the humidifier is on, this feature allows the airflow to continue for about 15 minutes at a low pressure (about 2 cmH ₂ O) after you press the Knob  to discontinue treatment. This will blow off the vapor left in the water chamber to avoid any damage to the device. When this feature is set to "Off," which means it is disabled, the airflow stops delivering air instantly after you press the Knob  .
Time	00:00 — 23:59	Set time by adjusting this option
Brightness	High / Low	Set screen brightness by adjusting this option.

Mask Type	Full Face; Nasal; Pillow; Other	There are three mask types available, Full Face (full-face mask), Nasal (nasal mask), and Pillow (nasal pillow mask). The patient can choose other suitable masks as well. When selecting masks other than the above three types of REACT HEALTH masks, the patient can identify the masks as other
iCode	iCode, iCode QR, iCode QR+	iCode provides access to the patient's compliance data during a recent time period. The iCode mode displays data in sequences of numbers, and the iCode QR / iCode QR+ mode displays data in two-dimensional codes
About	—	Displays related information of the device (Model, SN, Version, ID). This is read-only and cannot be edited

16. Alert

Alert Message	Description
Power Failure!!!	An audible alert will sound in 6s if the device is accidentally disconnected from power when it is delivering air. Note: (1) The alert will not sound if power failure occurs when the device is in standby state. (2) No alert message on the screen during a power failure
Device Fault!!!	An audible alert will sound if no airflow comes out of the machine; the screen will display " Device Fault!!! "
Leak!!	When the airflow is on, an audible alert will sound in 40 s if the air leak rate is excessive; the screen will display " Leak!! "
Low Input Voltage!!	If the voltage supplied by power adaptor is lower than 22V, an audible alert will sound and the screen will display " Low Input Voltage!! "
Humidifier Failure!!	When humidifier is applied, an audible alert will sound when the humidifier fails to work; the screen will display " Humidifier Failure!! "
Please Change Filter!	When the Filter Alert feature is enabled, an audible alert will sound if the preset replacement time reaches but without replacing the air filter; the screen will display " Please Change Filter! "
SD Card Full!	The screen will display " SD Card Full! " if the SD card has reached its maximum capacity
Reinsert SD card!	The screen will display " Reinsert SD card! " if the SD card fails to work

17. Cleaning

WARNINGS!

- Cleaning the device and its accessories as recommended in the following sections is essential to prevent respiratory infections.
- To avoid electric shock, always unplug the device before cleaning.
- Use mild soap that is nontoxic to humans.
- Follow the manufacturer's instructions on cleaning the mask and tube and on determining the frequency of cleaning.
- Before cleaning, check whether the device has been disconnected from the power supply, whether the power cord has been unplugged, and whether the water chamber of the device has cooled down. Make sure the plate has cooled down to room temperature, so you do not get burned.
- The device shall not be serviced or maintained while in use with a patient.
- Sterilization of this device and its components other than recommended is not permitted.
- Do not open or modify the device. There are no user serviceable parts inside. Repairs and servicing should only be performed by an authorized service agent.
- Avoid the use of any CPAP cleaner or disinfection device that relies on ozone (i.e. activated oxygen). Device warranty may terminate if the damage is caused by the use of an ozone cleaner.
- In order to prevent contamination of the device, use only manufacturer-approved filters on this device conforming to ISO 23328-1:2003 and ISO 23328-2:2002 standards.
- The device shall not be serviced or maintained while a patient is using it.

CAUTIONS!

- Overheating of the materials could lead to early fatigue of these materials.
- Do not use solutions containing chlorinated lime, chlorine, or aromatic to clean the device and its accessories. Liquid soap containing moisturizing agents or antimicrobials should not be used either. These solutions may harden cleaned materials or reduce their lifespan.
- Do not clean or dry the device and its accessories when the temperature is higher than 80°C (176°F). High temperatures could reduce product life.
- Do not immerse the device in any fluids.

17.1 Cleaning the Mask and Headgear

For details, refer to the cleaning instructions in the user manual for the mask.

17.2 Cleaning the Cellular Module

Wipe the surface of the Cellular Module with a clean, soft, and slightly damp cloth.

CAUTION!

- It is recommended to clean the Cellular Module once a week.

17.3 Cleaning the Water Chamber

Preparation tools: soft bristle brush, drinking quality water, mild liquid detergent.

Remove the Water Chamber according to instructions in 14.1.1.

(1) **Opening the Water Chamber:** Open the cap of the water chamber, as shown in Fig. 17-1.

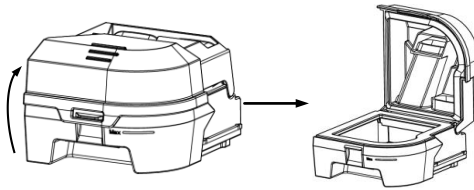


Fig. 17-1

(2) **Cleaning the Water Chamber:** Clean the water chamber with a soft bristled brush for one minute after soaking the water chamber in detergent solution (see table below) for 5 minutes. Pay particular attention to all crevices and cavities. Then rinse with running water for 5 minutes. Wipe it dry with a soft cloth or air dry out of direct sunlight.

Detergent	Water temperature
Alconox. (diluted at 1%)	Warm water (approx. 113 to 140°F or 45 to 60°C)

(3) **Returning the Water Chamber** according to instructions in 14.1.3.

WARNINGS!

- Emptying and cleaning the water chamber daily will help prevent mold and bacteria growth.
- Allow the water in the chamber to cool down to room temperature before removing it from the device.

CAUTIONS!

- Clean the water chamber only after the water in it cools. Make sure that no water enters the device.
- After cleaning, rinse the water chamber thoroughly in clean water to make sure that no soap residue is left; then wipe it dry with a lint-free cloth, so as to prevent calcareous accumulations.
- Check the water chamber for any leak or damage. Replace the water chamber if any damage is present.
- It is recommended to do daily cleaning of the water chamber.

17.4 Cleaning the Enclosure

Wipe the surface of the device with a soft, slightly damp cloth.

CAUTIONS!

- The device can only be used after the enclosure is dry, so that no moisture enters the device.
- It is recommended to clean the enclosure once a week.

17.5 Cleaning the Tube

For details, refer to the cleaning instructions in the user manual of the L1 tube.

CAUTION!

- It is recommended to clean the tube once a week.

17.6 Cleaning and Replacing the Air Filter

- (1) Open the air filter cap to remove the air filter.
- (2) Put the new air filter in the filter area, and then place the filter cap back properly.
- (3) When the filter is dirty, the filter can be cleaned as follows:
Hold the air filter and align it with the tap water flow. Rinse the air filter with running water. When cleaning, gently press the air filter, but do not pull the air filter.
Allow the filter to air dry completely before reinstalling it.

CAUTION: Never install a wet filter into the device. You must ensure sufficient drying time for the cleaned filter.

CAUTIONS!

- To avoid material damage, do not place the spare air filter in direct sunlight, humid environments, or temperatures below the freezing point. The air filter should be replaced every 6 months (It may be replaced more frequently based on actual sanitary conditions).
- Operating the device with a dirty air filter may stop it from working properly and may cause damage to the device.
- Replacing the air filter and filter cap, device must be unplugged.

18. Traveling

18.1 Traveling with the Device

(1) Use the REACT HEALTH carrying case to carry the device and accessories along with you. Do not put them in your checked baggage.

(2) This device operates on power supplies of 100 V - 240 V and 50 Hz / 60 Hz, and is suitable for use in any country in the world. No special adjustment is necessary, but you will need to find out the types of the power sockets in your destination. Bring, if necessary, a power socket adaptor which can be purchased in electronics stores.

(3) Remember to bring a spare air filter and the emergency documents (filled and signed by your physician) about this device. If you plan to travel by air, remember to bring the multi-language emergency documents about respiratory therapy, in case that the border and customs officers in your destination country inspect the device. With the emergency documents, you can prove to them that it is a medical device.

(4) Security Stations: For convenience at security stations, there is a note on the bottom of the device stating that it is medical equipment. It may be helpful to bring this manual along with you to help security personnel understand the device.

CAUTIONS!

- Empty the water chamber before packing the device for your trip; in order to prevent any remaining water from entering the device.
- If the device is used when the atmospheric pressure is out of the stated range (See Section 5), the accuracy of the leakage alert will be affected.

18.2 Traveling by airplane

For some airlines, medical devices do not count toward carry-on luggage limits. Please check with your airline for their policy regarding medical equipment.

You can use your device on a plane as it meets the Federal Aviation Administration (FAA) requirements.

Aircraft Use

REACT HEALTH confirms that the device meets the Federal Aviation Administration (FAA) requirements (RTCA DO 160, section 20, category T and section 21, category M) for all phases of air travel.

19. Reordering

Contact your home care provider to order accessories or replacement filters. The device does not require routine servicing.

WARNINGS!

- If you notice any unexplained changes in the performance of the device, if it is

making unusual or harsh sounds, if it has been dropped or mishandled, if the enclosure is broken, or if water has entered the enclosure, discontinue use. Contact your home care provider.

- If the device malfunctions, contact your home care provider immediately. Never attempt to open the enclosure of the device. Repairs and adjustments must be performed by REACT HEALTH -authorized service personnel only. Unauthorized service could cause injury, invalidate the warranty, or result in costly damage.
- If necessary, contact your local authorized dealer or REACT HEALTH for technical support and documents.

20. Technical Support

Please contact REACT HEALTH directly if you need the circuit diagram of the device and the list of components for certain purposes such as maintenance or connection to other equipment. REACT HEALTH will provide the circuit diagram and / or other technical documents in whole or in part according to your needs.

21. Disposal

When the device reaches the end of its service life, dispose of the device and packaging in accordance with local laws and regulations.

22. Troubleshooting

The table below lists common problems you may have with the device and possible solutions to those problems. If none of the corrective actions solve the problem, contact your home care provider.

22.1 Common Problems in Patients and Corresponding Solutions

Problem	Possible Cause	Solution (s)
Dry, cold, runny, and blocked nose; having a cold	The nose reacts to the airflow and cold. Due to fast airflow, the air becomes cold, leading to nasal mucosa irritation and subsequent dryness and swelling	Increase the humidity setting of the device. Contact your physician, and continue treatment unless the physician suggests otherwise
Dry mouth and throat	Probably because the patient sleeps with his or her mouth open, and the pressurized air goes out via the mouth, leading to nasal and throat dryness	Use a chin strap to prevent the mouth from opening during sleep, or use a full-face mask. Contact your physician for details
Eye irritation	The mask size or model may not be correct, or the mask is not positioned correctly, thereby leading to air leakage	Narrow the distance between the forehead support of the mask and the forehead. Note that adjusting the mask too tight may leave markings on the patient's face. Add additional filling to the mask so it does not leak. Contact your home care provider for an appropriate mask. Add additional filling to the mask if necessary
	Mask cushion (the soft part of the mask) hardens	Replace the mask or mask cushion
Facial reddening	The mask is too tight	Loosen the headgear
	The distance between the forehead support of the mask and the forehead is not correct	Try a different distance. The angle and size of the forehead support differ according to the type of masks

Problem	Possible Cause	Solution (s)
Facial reddening	Wrong mask size	Contract your home care provider for a correct-size mask
	The patient is allergic to the materials of the mask	Contact your physician and home care provider. Use a mask which is not made with natural rubber latex. Place a lining between the skin and mask
Water in mask	When the humidifier is used, the humidified air tends to condense in the cold tube and mask if the room temperature is low	Turn the humidity setting down, or raise the room temperature. Place the tube under the quilt, or use the tube cover. Hang the tube loosely, and the lowest part of the tube should be lower than the patient's head
Nasal, sinus, or ear pain	Sinus or middle ear inflammation	Contact your physician immediately
Discomfort due to inability to adapt to the treatment pressure	The patient will feel uncomfortable when the treatment pressure is higher than 13 cmH ₂ O. However, the treatment pressure is determined according to the patient's conditions, and cannot treat sleep apnea if the treatment pressure is set too low	It takes a maximum of four weeks to adapt to pressurized air. Relax and breathe through the nose. If the problem still exists, contact your physician
Obstructive sleep apnea symptoms recur	Probably because the patient sleeps with his or her mouth open, and the pressurized air goes out via the mouth, leading to blockage in the respiratory tract	Use a chin strap to prevent the mouth from opening during sleep, or use a full-face mask. Contact your physician for details
The device is too noisy	The tube is not connected properly	Reconnect the tube properly
Air delivered from the device is abnormally hot	The air inlet of the device may be partially blocked, leading to insufficient airflow into the device	Replace the air filter (see 17.6 Cleaning and Replacing the Air Filter), and clean the air inlet
		Place the device in an area where air flows freely, and make sure the device is at least 20 centimeters away from the wall, curtain, or other things


22.2 Common Problems in the Device and Corresponding Solutions

Problem	Possible Cause	Solution (s)
The device does not work when it is turned on	The Auto On / Off feature is enabled	Take a few deep breaths with the mask on, and the device will start automatically
	Power is not connected properly	Ensure that the power cord, power adapter, and the device are connected properly
	There is no voltage	Check whether a power outage occurs by turning on a light or other means. If you are sure the fuse in the device is broken, contact your home care provider for repair
	Cannot find any cause	Contact your home care provider
The device is working, but the pressure inside the mask differs from the set treatment pressure	The tube is not connected properly	Reconnect the tube properly
	There may be holes in the mask or pressure sensing tube	Contact your home care provider
	It is a faulty device	Contact your home care provider
The device produces very low pressures	The air inlet of the device may be blocked	Replace the air filter (see 17.6 Cleaning and Replacing the Air Filter), and clean the air inlet. Make sure the air inlet is unblocked
	The treatment pressure has been changed accidentally	Contact your physician
	When the Ramp feature is enabled, it takes some time for the initial pressure to rise to the treatment pressure. This is normal	If necessary, disable the Ramp feature, or set the ramp time shorter
After the device is turned on, the screen displays intermittently, or displays nothing at all	The operating system of the device needs to be readjusted or restarted	Unplug the power cord of the device, and re-plug it 20 seconds later
The device is in standby, and will not start	The operating system of the device needs to be readjusted or restarted	Unplug the power cord of the device, and re-plug it 20 seconds later

23. EMC Requirements

Guidance and manufacturer's declaration - electromagnetic emissions		
The device is intended for use in the electromagnetic environment specified below. The user of the device should ensure that it is used in such an environment.		
Emissions Test	Compliance	Electromagnetic Environment - Guidance
RF emissions CISPR 11	Group 1	The device uses RF energy only for its internal function. Therefore its RF emissions are very low and are not likely to cause any interference in nearby electronic equipment
RF emissions CISPR 11	Class B	The device is suitable for use in all establishments including domestic establishments and those directly connected to the public low-voltage power supply network that supplies buildings used for domestic purposes
Harmonic emissions IEC 61000-3-2	Class A	
Voltage fluctuations / flicker emissions IEC 61000-3-3	Complies	

Guidance and manufacturer's declaration - electromagnetic immunity			
The device is intended for use in the electromagnetic environment specified below. The user of the device should make sure that it is used in such an environment.			
Immunity Test	IEC 60601 Test Level	Compliance Level	Electromagnetic Environment - Guidance
Electrostatic discharge (ESD) IEC 61000-4-2	±8 kV contact ±15 kV air	±8 kV contact ±15 kV air	Floor should be wood, concrete or ceramic tile. If floors are covered with synthetic material, the relative humidity should be at least 30%
Electrical fast transient / burst IEC 61000-4-4	±2 kV for power supply lines	±2 kV for power supply lines	Mains power quality should be that of a typical commercial or hospital environment
Surge IEC 61000-4-5	±1 kV Line (s) to line (s)	±1 kV Line (s) to line (s)	Mains power quality should be that of a typical commercial or hospital environment
Voltage dips, short interruptions and voltage variations on power supply input lines IEC 61000-4-11	0% U_T ; 0.5 cycle At 0°, 45°, 90°, 135°, 180°, 225°, 270° and 315° 0% U_T ; 1 cycle 70% U_T ; 25 / 30 cycle At 0° 0% U_T ; 250 / 300 cycle	0% U_T ; 0.5 cycle At 0°, 45°, 90°, 135°, 180°, 225°, 270° and 315° 0% U_T ; 1 cycle 70% U_T ; 25 / 30 cycle At 0° 0% U_T ; 250 / 300 cycle	Mains power quality should be that of a typical commercial or hospital environment. If the user of the device requires continued operation during power mains interruptions, it is recommended that the device be powered from an uninterruptible power supply or a battery
Power frequency (50 Hz / 60 Hz) magnetic field IEC 61000-4-8	30 A/m	30 A/m	Power frequency magnetic fields should be at levels characteristic of a typical location in a typical commercial or hospital environment
Note: U_T is the AC mains voltage prior to application of the test level.			

Guidance and manufacturer's declaration - electromagnetic immunity			
The device is intended for use in the electromagnetic environment specified below. The user of the device should make sure that it is used in such an environment.			
Immunity Test	IEC 60601 Test Level	Compliance Level	Electromagnetic Environment - Guidance
Conducted RF IEC 61000-4-6	3 V 0.15 MHz ~ 80 MHz 6 V in ISM and amateur radio bands between 0.15 MHz and 80 MHz	3 V 0.15 MHz ~ 80 MHz 6 V in ISM and amateur radio bands between 0.15 MHz and 80 MHz	Portable and mobile RF communications equipment should be used no closer to any part of the device, including cables, than the recommended separation distance calculated from the equation applicable to the frequency of the transmitter. Recommended separation distance $d = 1.17\sqrt{P}$ $d = 0.35\sqrt{P}$ 80 MHz to 800 MHz $d = 0.70\sqrt{P}$ 800 MHz to 2.5 GHz Where P is the maximum output power rating of the transmitter in watts (W) according to the transmitter manufacturer and d is the recommended separation distance in meters (m). Field strengths from fixed RF transmitter, as determined by an electromagnetic site survey, ^a should be less than the compliance level in each frequency range. ^b Interference may occur in the vicinity of equipment marked with the following symbol: 
Radiated RF IEC 61000-4-3	10 V/m 80 MHz to 2.7 GHz	10 V/m 80 MHz to 2.7 GHz	
Note 1: At 80 MHz and 800 MHz, the higher frequency range applied.			
Note 2: These guidelines may not apply in all situations. Electromagnetic propagation is affected by absorption and reflection from structures, objects and people.			
^a Field strengths from fixed transmitters, such as base stations for radio (cellular / cordless) telephones and land mobile radios, amateur radio, AM and FM radio broadcast and TV broadcast cannot be predicted theoretically with accuracy. To assess the electromagnetic environment due to fixed RF transmitters, an electromagnetic site survey should be considered. If the measured field strength in the location in which the device is used exceeds the applicable RF compliance level above, the device should be observed to verify normal operation. If abnormal performance is observed, additional measures may be necessary, such as re-orienting or relocating the device.			
^b Over the frequency range 150 kHz to 80 MHz, the field strengths should be less than 10 V/m.			

Recommended separation distances between portable and mobile RF communications equipment and the device

The device is intended for use in an electromagnetic environment in which radiated RF disturbances are controlled. The customer or the user of the device can help prevent electromagnetic interference by maintaining a minimum distance between portable and mobile RF communications equipment (transmitters) and the device as recommended below, according to the maximum output power of the communications equipment.

Rated maximum output of transmitter (W)	150 kHz ~ 80 MHz $d = 1.17\sqrt{P}$	80 MHz ~ 800 MHz $d = 0.35\sqrt{P}$	800 MHz ~ 2.5 GHz $d = 0.70\sqrt{P}$
0.01	0.12	0.04	0.07
0.1	0.37	0.12	0.23
1	1.17	0.35	0.70
10	3.70	1.11	2.22
100	11.7	3.50	7.00


Note 1: At 80 MHz and 800 MHz, the higher frequency range applied.

Note 2: These guidelines may not apply in all situations. Electromagnetic propagation is affected by absorption and reflection from structures, objects and people.

For transmitters rated at a maximum output power not listed above, the recommended separation distance d in meters (m) can be estimated using the equation applicable to the frequency of the transmitter, where P is the maximum output power rating of the transmitter in watts (W) according to the transmitter manufacturer.

Recommended separation distances between RF wireless communications equipment

The device is intended for use in an electromagnetic environment in which radiated RF disturbances are controlled. The customer or the user of the device can help prevent electromagnetic interference by maintaining a minimum distance between RF wireless communications equipment and the device as recommended below, according to the maximum output power of the communications equipment.

Frequency (MHz)	Maximum Power (W)	Distance	IEC 60601 Test Level	Compliance Level	Electromagnetic Environment - Guidance
385	1.8	0.3	27	27	RF wireless communications equipment should be used no closer to any part of the device, including cables, than the recommended separation distance calculated from the equation applicable to the frequency of the transmitter. Recommended separation distance $E = \frac{6}{d} \sqrt{P}$ Where P is the maximum output power rating of the transmitter in watts (W) according to the transmitter manufacturer and d is the recommended separation distance in meters (m). Field strengths from fixed RF transmitter, as determined by an electromagnetic site survey, should be less than the compliance level in each frequency range. Interference may occur in the vicinity of equipment marked with the following symbol: 
450	2	0.3	28	28	
710	0.2	0.3	9	9	
745					
780					
810					
870	2	0.3	28	28	
930	2	0.3	28	28	
1720					
1845					
1970					
2450					
5240	0.2	0.3	9	9	
5500					
5785					

Note: These guidelines may not apply in all situations. Electromagnetic propagation is affected by absorption and reflection from structures, objects and people.

WARNINGS!

- The device should not be used in the vicinity or on the top of other electronic equipment such as cell phone, transceiver or radio control products. If you have to do so, the device should be observed to verify normal operation.
- Use of this equipment adjacent to or stacked with other equipment should be avoided because it could result in improper operation. If such use is necessary, this equipment and the other equipment should be observed to verify that they are operating normally.
- Use of accessories, transducers and cables other than those specified or provided by the manufacturer of this equipment could result in increased electromagnetic emissions or decreased electromagnetic immunity of this equipment and result in improper operation.
- Portable RF communications equipment (including peripherals such as antenna cables and external antennas) should be used no closer than 30 cm (12 inches) to any part of the Luna® II CPAP / Luna® II Auto CPAP, including cables specified by the manufacturer. Otherwise, degradation of the performance of this equipment could result.
- The device may be interfered with by other equipment, even if that other equipment complies with CISPR EMISSION requirements.
- The device may be subject to interference by the electromagnetic field of some known or unknown radio frequency transmitters in the environment during use. If interference occurs, please stay away from the interfered electromagnetic environment, or find and turn off the electromagnetic field interference source before continuing to use it.
- When the device is exposed to welding, electrosurgery, defibrillation, X-ray (γ ray), infrared radiation, transient electromagnetic field, including nuclear magnetic resonance (MRI) and radio interference environment, the product may be damaged.
- During operation of the device, due to electrostatic interference, the following phenomena may occur: (1) Temporary loss of function or degradation of performance, such as abnormal screen display, etc. The device will recover to normal after being restarted; (2) Automatic restart of the device. These phenomena will not affect the normal use of the device, and will not cause permanent performance degradation or function loss of the device.

24. Limited Warranty

REACT HEALTH warrants that the device shall be free from defects of workmanship and materials and will perform in accordance with the product specifications for a period of two (2) years for main unit and three (3) months for all accessories from the date of sale by REACT HEALTH to the dealer.

If the product fails to perform in accordance with the product specifications, REACT HEALTH will repair or replace, at its option, the defective material or part. REACT HEALTH will pay customary freight charges from REACT HEALTH to the dealer location only. This warranty does not cover damage caused by accident, misuse, abuse, alteration and other defects not related to material or workmanship.

REACT HEALTH DISCLAIMS ALL LIABILITY FOR ECONOMIC LOSS, LOSS OF PROFITS, OVERHEAD OR CONSEQUENTIAL DAMAGES WHICH MAY BE CLAIMED TO ARISE FROM ANY SALE OR USE OF THIS PRODUCT. SOME STATES DO NOT ALLOW THE EXCLUSION OR LIMITATION OF INCIDENTAL OR CONSEQUENTIAL DAMAGES, SO THE ABOVE LIMITATION OR EXCLUSION MAY NOT APPLY TO YOU.

To exercise the rights under this warranty, contact the local authorized dealers or:

REACT HEALTH

5101 Fruitville Rd., Suite 200
Sarasota, FL 34232
T: (863) 226-6285

For additional information, please visit our Patient Portal at:

www.reacthealth.com
icodeconnect.com – Web-based cloud for report generation and storage

BMC Medical Co., Ltd.

Room 10, 17F, Building 4, Huiya Plaza, No.16 Lize Road, Fengtai District, 100073 Beijing,
PEOPLE'S REPUBLIC OF CHINA
Tel: +86-10-51663880
Fax: +86-10-51663880 Ext. 810
URL: en.bmc-medical.com
E-mail: intl@bmc-medical.com